



Human zinc transporter 8 (SLC30A8) ELISA Kit

Catalog Number. CSB-EL021573HU

For the quantitative determination of human zinc transporter 8 (SLC30A8) concentrations in serum, plasma, tissue homogenates, cell lysates.

This package insert must be read in its entirety before using this product.

If You Have Problems

Technical Service Contact information

Phone: 86-27-87582341

Fax: 86-27-87196150

Email: tech@cusabio.com

Web: www.cusabio.com

In order to obtain higher efficiency service, please ready to supply the lot number of the kit to us (found on the outside of the box).

PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. Antibody specific for SLC30A8 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells with a Horseradish Peroxidase (HRP) conjugated antibody specific for SLC30A8. Following a wash to remove any unbound reagent, a substrate solution is added to the wells and color develops in proportion to the amount of SLC30A8 bound in the initial step. The color development is stopped and the intensity of the color is measured.

DETECTION RANGE

80 pg/ml-1400 pg/ml.

SENSITIVITY

The minimum detectable dose of human SLC30A8 is typically less than 50 pg/ml.

The sensitivity of this assay, or Lower Limit of Detection (LLD) was defined as the lowest human SLC30A8 concentration that could be differentiated from zero. It was determined the mean O.D value of 20 replicates of the zero standard added by their three standard deviations.

SPECIFICITY

This assay has high sensitivity and excellent specificity for detection of human SLC30A8. No significant cross-reactivity or interference between human SLC30A8 and analogues was observed.

Note: Limited by current skills and knowledge, it is impossible for us to complete the cross-reactivity detection between human SLC30A8 and all the analogues, therefore, cross reaction may still exist.

PRECISION

Intra-assay Precision (Precision within an assay): CV%<15%

Three samples of known concentration were tested twenty times on one plate to assess.

Inter-assay Precision (Precision between assays): CV%<15%

Three samples of known concentration were tested in twenty assays to assess.

LIMITATIONS OF THE PROCEDURE

- **FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.**
- The kit should not be used beyond the expiration date on the kit label.
- Do not mix or substitute reagents with those from other lots or sources.
- If samples generate values higher than the highest standard, dilute the samples and repeat the assay.
- Any variation in operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
- This assay is designed to eliminate interference by soluble receptors, binding proteins, and other factors present in biological samples. Until all factors have been tested in the Immunoassay, the possibility of interference cannot be excluded.

MATERIALS PROVIDED

Reagents	Quantity
Assay plate	1 (96 wells)
Standard	6 (lyophilized)
HRP-conjugate	1 x 6 ml
Wash Buffer (20 x concentrate)	1 x 15 ml
Substrate A	1 x 7 ml
Substrate B	1 x 7 ml
Stop Solution	1 x 7 ml
Adhesive Strip (For 96 wells)	4
Instruction manual	1

STANDARD CONCENTRATION

Standard	S0	S1	S2	S3	S4	S5
Concentration (pg/ml)	0	80	160	320	800	1400

STORAGE

Unopened kit	Store at 2 - 8°C. Do not use the kit beyond the expiration date.
Opened kit	May be stored for up to one month at 2 - 8° C.

***Provided this is within the expiration date of the kit.**

OTHER SUPPLIES REQUIRED

- Microplate reader capable of measuring absorbance at 450 nm, with the correction wavelength set at 600 nm - 630 nm.
- An incubator which can provide stable incubation conditions up to 37°C±0.5°C.
- Squirt bottle, manifold dispenser, or automated microplate washer.
- Absorbent paper for blotting the microtiter plate.
- 100 mL and 500 mL graduated cylinders.
- Deionized or distilled water.

- Pipettes and pipette tips.
- Test tubes for dilution.

PRECAUTIONS

The Stop Solution provided with this kit is an acid solution. Wear eye, hand, face, and clothing protection when using this material.

SAMPLE COLLECTION AND STORAGE

- **Serum** Use a serum separator tube (SST) and allow samples to clot for two hours at room temperature or overnight at 4°C before centrifugation for 15 minutes at 1000 xg. Remove serum and assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.
- **Plasma** Collect plasma using EDTA, or heparin as an anticoagulant. Centrifuge for 15 minutes at 1000 xg at 2-8°C within 30 minutes of collection. Assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.
- **Tissue Homogenates** 100mg tissue was rinsed with 1X PBS, homogenized in 1 ml of 1X PBS and stored overnight at -20°C. After two freeze-thaw cycles were performed to break the cell membranes, the homogenates were centrifuged for 5 minutes at 5000 x g, 2 - 8°C. The supernate was removed and assayed immediately. Alternatively, aliquot and store samples at -20°C or -80°C. Centrifuge the sample again after thawing before the assay. Avoid repeated freeze-thaw cycles.
- **Cell Lysates (1) Adherent Cell:** Remove media and rinse cells once with ice-cold PBS (pH7.2-7.4). Scrape cells off the plate and transfer to an appropriate tube. Dilute cell suspension with 1xPBS (pH7.2-7.4), until cell concentration reached 100 million/ml. Then store overnight at -20°C. After two freeze-thaw cycles to break up the cell membranes, the cell lysates were centrifuged for 5 minutes at 5000 x g, 2 - 8°C. Collect the

supernatant. Cell lysates should be assayed immediately or aliquotted and stored at -20°C. Centrifuge the sample again after thawing before the assay. Avoid repeated freeze-thaw cycles. **(2) Suspension Cell:** Collect cells with appropriate tube, centrifuge for 5 minutes at 1000x g, 2 - 8°C. Remove the supernatant and resuspend cells with 1xPBS (pH7.2-7.4). Centrifuge for 5 minutes at 1000 x g, 2 - 8°C. Remove the supernatant. Dilute cell with 1xPBS (pH7.2-7.4), until cell concentration reached 100 million/ml. Then store overnight at -20°C. After two freeze-thaw cycles to break up the cell membranes, the cell lysates were centrifuged for 5 minutes at 5000 x g, 2 - 8°C. Collect the supernatant. Cell lysates should be assayed immediately or aliquotted and stored at -20°C. Centrifuge the sample again after thawing before the assay. Avoid repeated freeze-thaw cycles.

Note:

1. CUSABIO is only responsible for the kit itself, but not for the samples consumed during the assay. The user should calculate the possible amount of the samples used in the whole test. Please reserve sufficient samples in advance.
2. Samples to be used within 5 days may be stored at 2-8°C, otherwise samples must be stored at -20°C (≤ 1 month) or -80°C (≤ 2 month) to avoid loss of bioactivity and contamination.
3. Grossly hemolyzed samples are not suitable for use in this assay.
4. If the samples are not indicated in the manual, a preliminary experiment to determine the validity of the kit is necessary.
5. Please predict the concentration before assaying. If values for these are not within the range of the standard curve, users must determine the optimal sample dilutions for their particular experiments.
6. Tissue or cell extraction samples prepared by chemical lysis buffer may cause unexpected ELISA results due to the impacts of certain chemicals.
7. Owing to the possibility of mismatching between antigen from other resource and antibody used in our kits (e.g., antibody targets conformational epitope rather than linear epitope), some native or recombinant proteins from other manufacturers may not be recognized by our products.
8. Influenced by the factors including cell viability, cell number and also sampling time, samples from cell culture supernatant may not be detected by the kit.
9. Fresh samples without long time storage are recommended for the test. Otherwise, protein degradation and denaturalization may occur in those samples and finally lead to wrong results.

REAGENT PREPARATION

Note:

- **Kindly use graduated containers to prepare the reagent.**
 - Bring all reagents to room temperature (18-25°C) before use for 30min.
 - Distilled water is recommended to be used to make the preparation for reagents. Contaminated water or container for reagent preparation will influence the detection result.
- 1、 **Wash Buffer(1x)**- If crystals have formed in the concentrate, warm up to room temperature and mix gently until the crystals have completely dissolved. Dilute 15 ml of Wash Buffer Concentrate (20 x) into deionized or distilled water to prepare 300 ml of Wash Buffer (1 x).

2、 **Standard**

Centrifuge the standard vial at 6000-10000rpm for 30s.

Reconstitute each lyophilized **Standard(S0-S5)** with **0.5 ml** of **ddH₂O**. Mix the standard to ensure complete reconstitution and allow the standard to sit for a minimum of 10 minutes with gentle agitation prior to use.

ASSAY PROCEDURE


Bring all reagents and samples to room temperature before use. Centrifuge the sample again after thawing before the assay. It is recommended that all samples and standards be assayed in duplicate.

1. Prepare all reagents, working standards, and samples as directed in the previous sections.
2. Determine the number of wells to be used and put any remaining wells and the desiccant back into the pouch and seal the ziploc, store unused wells at 4°C.
3. Set a **Blank** well without any solution.
4. Add 50µl of **Standard** or **Sample** per well. Standard need test in duplicate.
5. Add 50µl of **HRP-conjugate** to each well(Not to **Blank** well). Mix well and then incubate for one hour at 37°C.
6. Aspirate each well and wash, repeating the process two times for a total of three washes. Wash by filling each well with **Wash Buffer** (200µl) using a squirt bottle, multi-channel pipette, manifold dispenser, or autowasher, and let it stand for 10 seconds, complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
7. Add 50µl of **Substrate A** and 50µl of **Substrate B** to each well, mix well. Incubate for 15 minutes at 37°C. Keeping the plate away from drafts and other temperature fluctuations in the dark.
8. Add 50µl of **Stop Solution** to each well, gently tap the plate to ensure thorough mixing.
9. Determine the optical density of each well within 10 minutes, using a microplate reader set to 450 nm.

Note:

1. The final experimental results will be closely related to validity of the products, operation skills of the end users and the experimental environments.
2. Samples or reagents addition: Please carefully add samples to wells and mix gently to avoid foaming. Do not touch the well wall as possible. For each step in the procedure, total dispensing time for addition of reagents or samples to the assay plate should not exceed 10 minutes. This will ensure equal elapsed time for each pipetting step, without interruption. Duplication of all standards and specimens, although not required, is recommended. To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
3. Incubation: To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary. Do not allow wells to sit uncovered for extended periods between incubation steps. Once reagents have been added to the well strips, DO NOT let the strips DRY at any time during the assay. Incubation time and temperature must be observed.
4. Washing: The wash procedure is critical. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading. When using an automated plate washer, adding a 30 second soak period following the addition of wash buffer, and/or rotating the plate 180 degrees between wash steps may improve assay precision.
5. Controlling of reaction time: Observe the change of color after adding Substrates (e.g. observation once every 10 minutes). Substrates should change from colorless or light blue to gradations of blue. If the color is too deep, add Stop Solution in advance to avoid excessively strong reaction which will result in inaccurate absorbance reading.
6. Substrates are easily contaminated. Substrates should remain colorless or light blue until added to the plate. Please protect it from light.
7. Stop Solution should be added to the plate in the same order as the Substrates. The color developed in the wells will turn from blue to yellow upon addition of the Stop Solution. Wells that are green in color indicate that the Stop Solution has not mixed thoroughly with the Substrates.

ASSAY PROCEDURE SUMMARY

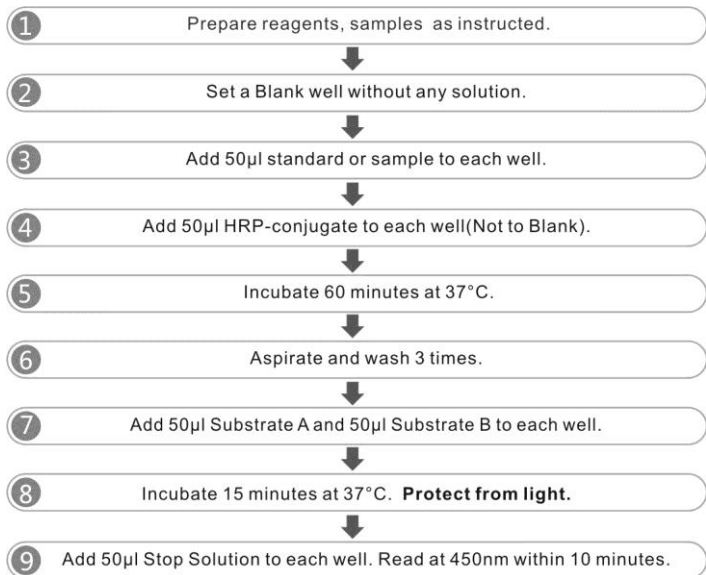

CALCULATION OF RESULTS

Using the professional soft "Curve Expert" to make a standard curve is recommended, which can be downloaded from our web.

Average the duplicate readings for each standard and sample and subtract the average optical density of Blank.

Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the x-axis against the concentration on the y-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the SLC30A8 concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.